

Gutter Crutch | Technical Spec Sheet

Medical Crutches | Aluminium Gutter Crutch


Product Overview

Description:	Gutter crutch, comfortable for medium to long term use
Type:	Forearm, adjustable
Tube finish:	Anodised matt 15 microns (μ)
Colours:	Available in natural anodised finish only
Plastics:	Black
Structure:	Adjustable aluminium tube / Steel forearm support gutter plate with sponge padding and velcro strap.
Leg adjustment:	Adjustable to 10 different lengths (25mm apart)
Forearm adjustment:	Adjustable to 5 different lengths (25mm apart)
Punch options:	Single; or double by arrangement

Classification: ISO 9999
Code: 120306

Key

1. Gutter Top Support Line
 2. Gutter Top
 3. Foot
 4. Leg Section Axis
 5. Handle
 6. Gutter Support Line
 7. Front Handle Reference Point
 8. Rear Handle Reference Point
- a . Arm Section Length
 h . Handle Length
 l . Leg Section Length
 α . Support Angle
 β . Grip Angle

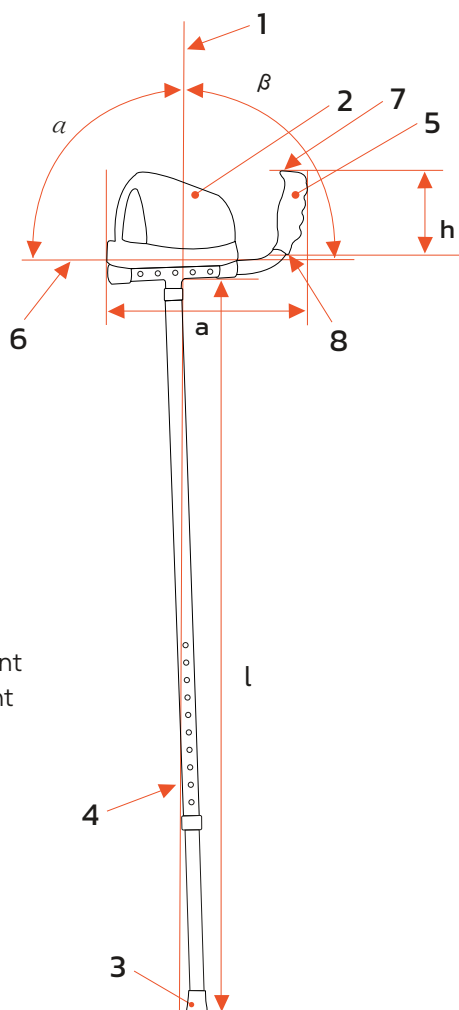

General Information

Support Angle (α)	90°
Grip Angle (β)	90°
Arm Section Length (a)	100 - 200mm
Mass	1145g
Max. Weight of User	100kg (pair)
Leg Section Length (l)	1100- 1200mm

Fittings

Cuff	Padded with sponge and velcro strap
Handle	Made of polypropylene Moulded handle (0.380 x 21.6mm)
Foot	Design for maximum adherence Fitted insidewith anti-perforation ring Replaceable

Support (Aluminium 6063-T6)

Length	750mm
Outside Diameter	25.4mm
Mass	253g
Height Adjustment	X10 holes
Plastics	Fitted with plastic slide at lower end for noise and vibration reduction

Lower Tube (Aluminium 6063-T6)

Foot Length	550mm
Outside Diameter	21.6mm
Mass	178g
Height Adjustment	Fitted with plated spring clip
Plastics	Fitted with plug at upper end for noise and vibration reduction
Foot	Fitted with 22mm foot plug at lower end

Label Details: Batch No. and Production Date; Maximum Weight = 100kg; Name and Address of Manufacturer
Tolerances: Allow for approximately 2% deviation of all stated measurements